



Pharmaworks, Inc., a Pro-Mach Company, following the guidelines under the Emergency Use Authorization (EUA) letter issued by the U.S. Food and Drug Administration on April 9, 2020, is now manufacturing face shields for health care professionals and first responders in accordance with CDC recommendations. Fully covering the front and sides of the face, this face shield does not impede head movement, allowing wearers a complete range of motion while providing a barrier of protection against droplets and splashes during the COVID-19 pandemic.

Designed for Flexibility

- Features a reusable, versatile rigid frame design.
- Die-cut Face Shield is ergonomically contoured.
- Face Shield can be cleaned and/or exchanged.
- A standard 3 hole punch can be used to make temporary replacement shields.
- No straps – uses a simple replaceable antimicrobial rubber band.
- Lightweight and comfortable to wear.

Specifications

- Molded Face Shield Frame: Polybutylene Thermoplastic, alcohol and solvent resistant
- Laser Cut Face Shield Frame: FDA Compliant Acrylic
- Face Shield Die Cut Sheet: PVC Film
- Face Shield Frame Cushion: Silicone or SBR Rubber
- Shield Dimensions: 12" wide x 8.5" tall for enhanced coverage.
- ANSI Z87.1-2015 Certified
- Single User Device, Reusable After Cleaning


Cleaning Recommendations

Frame

The face shield frame can be wiped down with alcohol or solvent-based cleaning solutions.

Face Shield

The face shield can be cleaned using soapy water and micro-fiber cloth. Avoid the use of ammonia based products as they contain harmful chemicals that can damage the shield's surface leaving it cloudy.

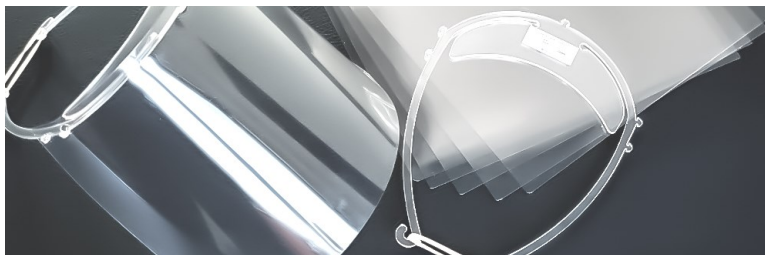

This Face Shield has not been approved by the FDA. The Face Shield has been authorized by the FDA under the EUA for use by healthcare providers as personal protective equipment. This Face Shield provides the minimal protection to the user during the COVID-19 pandemic when used consistently. This product is only authorized for the duration of the declaration that circumstances justifying the authorization of emergency use under Section 564(b)(1) of the Act, 21 USC 360bbb-3(b)(1) unless the authorization is terminated or revoked sooner. Face Shield should not be used in the presence of high intensity heat source or flammable gas.

Die Cut, Ergonomically Contoured Face Shield



Reusable Molded Frame

